**EN**

With your inspiring rating, COSTWAY will be more consistent to offer you EASY SHOPPING EXPERIENCE, GOOD PRODUCTS and EFFICIENT SERVICE!

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First, we appreciate your support. Please read the instructions carefully before operating this machine, and well keep it for later use. Note: we reserve the right to revise or explain the design of this product without further notice.

## FOOT MASSAGER SPECIFICATIONS

○ MODLE : LF-23

○ VOLTAGE : AC100~240V/50Hz 1,6A

○ POWER : 80W

○ RATED WORKING TIME : 15 MINUTES

## A/C ADAPTOR

○ INPUT : AC100~240V/50Hz 1,6A


## FEATURES

1. An entirely new massager specifically designed to massage your feet, ankles and calves all at the same time.
2. 4 parts of kneading discs in one appliance to massage you feet, ankles and calves.
3. Flexible kneading discs create a professional massage experience.
4. Strong kneading action penetrates deeply into the muscles of you feet, ankles and calves.
5. Pressure nodes located under your feet stimulate the reflex points on the soles of your feet.
6. Personal-preference of 3 kneading massage programs (P1/P2/P3)
7. The synergistic combination of kneading and vibratory action gives you a most relaxing and revitalizing massage experience, while achieving a beneficial reflexology and beautifying effect.
8. Strategically located buttons on the switch panel allow you to control the massage function easily.
9. Auto shut-off feature activates after the massager has been used for 15 minutes.
10. A pair of detachable fabric cover(with zipper)allows for cleaning to maintain good hygiene.
11. COULD accelerate circulation after USING THE HEATED FUNCTION ACCORDING TO YOUR NEED.

## EFFECTIVENESS


- 1.The kneading and vibratory actions go deep into the different layers of muscles of your feet, ankles, and calves, allowing you to rejuvenate and relax.
- 2.Kneading and vibration massage target vital reflex points on the feet and soles, achieving positive reflexology effects.
- 3.Stimulates blood circulation, thereby preparing and warming the muscles.
- 4.Relieves tension and reduces stress.
- 5.Relieves sore and tired feet effectively.
- 6.Helps to tone, shape and beautify your legs.

## NAMES OF PARTS AND FUNCTIONS

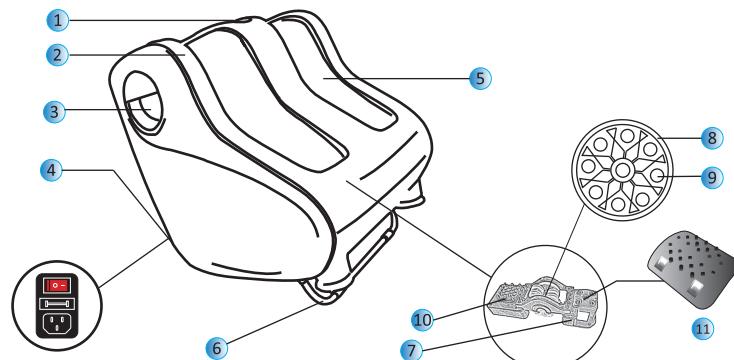- 1.Computer panel
- 2.Shell of mainframe
- 3.Handle
- 4.Power switch ( Back )
- 5.Cover ( with zipper )
- 6.Bracket
- 7.Roller assembly
- 8.Roll ball
- 9.Magnet
- 10.Heat up part
- 11.Vibration motor

## Operation Instruction of computer Panel

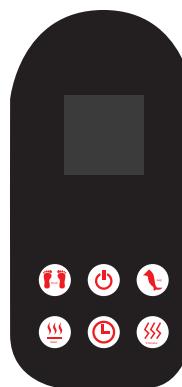
- |  |              |
|--|--------------|
|  | Power button |
|  | Foot massage |
|  | Calf massage |
|  | Time control |
|  | Vibrate      |
|  | Heat         |

### To activate the kneading discs only

- 1.Place the lower section of your legs into the fabric-covered slots and press ) to start.
- 2.Press the button and the default mode will be " ", press it again or twice, and mode will switch to " " and " ".
- 3.Press the button and the default mode will be " ", press it again or twice, and mode will switch to " " and " ".
- 4.Press the button and the default program will be "Vibrate", and press again to close.
- 5.Press the button and the default program will be "Heat", press again to close.
- 6.Press the button again to stop all of the programs.
- 7.Press the button to adjust the time of using.

Press the power switch in the back of the massager, the indicating light of stand-by will turn on.

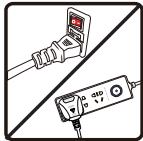Insert the adapter into the back seat of the massager in the correct direction. After plugging in the power source, the adapter indicator will be on.


Open or close bracket according to the height of feet.

Choose a chair that matches your height and place your feet in the massage trough.

Stop use: please press the off button, turn off the power switch, then cut the power and remove the adapter.

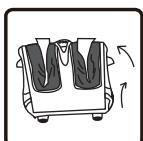If the power is interrupted during use, the first step is to put your feet on their side and move them easily.

If the power is interrupted during use, the second step: raise your knees.

## Caution

- Do not stand on, or exert excessive pressure on the massager. This may significantly slows down or stops the motor, potentially causes overheating and subsequent damage to the massager.
- The total recommended usage time of the massager is 15 minutes.
- Press the "stop" button to terminate massager if discomfort is experienced. The 4 pairs of discs will not stop instantaneously; they will continue to complete the last kneading cycle. Before stopping at their "home" position.
- Some users may experience itch on their feet during and after the massage. This is alright because it is the result of improved blood circulation.

## SAFETY

### Usage Environment

Do not use the massager in an excessively humid environment such as the bathroom.

Do not use the massager on a wet body, this may result in a short circuit or an electric shock.

### Suitability of User

The following users must consult their doctors before use.

Those suffering from cancer.

Those suffering from osteoporosis.

Those suffering from malignant tumors.

Those who have a high fever.

Women who are pregnant or menstruating.

### Safety Precautions during Use

Do not use the massager for any purpose, other than massage, this may result in accident or injury.

Do not use the massager when the fabric cover is torn, or parts of it are missing.

Do not let children, or persons who are not capable of handling equipment to operate the massager. Please keep the massager out of reach of children.

Do not stand or sit on the massager.

Do not insert metallic objects, or flammable materials into the main body. This may result in electrocution, short circuit and possible fire.

Do not use the massager within an hour of mealtime.

Do not cover the massager during operation. This may result in a rise in motor temperature and subsequent breakdown or fire.

Do not unplug switch off the massager during operation. This may result in injury.

Do not sleep while using the massager.

#### **Instances when the massager should not be used**

In the following events, please switch off the power immediately and disconnect the AC power cord from the electric outlet.

If water is accidentally spilled on the massager, this may result in electric shock.

If severe or abnormal pain is felt during massage, stop operation immediately and consult your doctor.

If you detect any malfunction or any other abnormal condition during operation.

If there is a power failure, injury may occur when power is restored unexpectedly.

If there is lightning.

#### **Assemble and repair the massager**

The massager must only be serviced or repaired by the service centre. Users must not disassemble or repair the massager themselves.

#### **Notes about power plug and cord**

Please check whether the voltage of power supply meets the adapter standard.

Do not disconnect the AC plug from the electrical outlet roughly.

Do not connect or disconnect the AC plug from the electrical outlet with wet hands.

It may result in electric shock or machine malfunction.

When disconnecting the AC power plug, pull it out by holding the plug instead of the power cord.

Do not place the AC cord under the massager or any heavy object.

Do not wind the AC cord around the massager as it may damage the cord and result in fire or electric shock.

Do not operate the massager with a damaged cord or tamper with the cord, contact the servicecentre for inspection.

Please use the adapter and attachment we recommend.

Do not use the massager if the electrical outlet is excessively loose.

#### **Product care and storage**

Be sure to remove the AC plug from the electrical outlet before cleaning the massager.

Do not use benzene, thinner or any other abrasive solvent on your massager.

Use a dry cloth when cleaning the massager. Never use a damp cloth to clean the massager or the power switch.

When necessary, detach the fabric covers for cleaning to maintain good hygiene.

When not in use for a long period of time, it is advisable that the AC power cord of the massager be coiled and the massager stored in a dust and moisture free environment.

Do not store the massager near intense heat or an open flame. Do not leave the massager under direct sunlight for long period of time.

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