



COSTWAY®

**Please give us a chance to make it right and do better !**

Contact our friendly customer service department for help first.  
Replacements for missing or damaged parts will be shipped ASAP !

**US** office: Fontana    **UK** office: Ipswich    **AU** office: Truganina  
**DE** office: Hamburg    **FR** office: Saint Vigor d'Ymonville    **PL** office: Gdańsk



## Contact Us!

**Do NOT return this item.**

Contact our friendly customer service department for help first.



US:cs.us@costway.com  
UK:cs.uk@costway.com  
AU:cs.au@costway.com

## USER'S MANUAL

### Rolling Bedside Commode JH10006

THIS INSTRUCTION BOOKLET CONTAINS **IMPORTANT** SAFETY INFORMATION. PLEASE READ AND KEEP FOR FUTURE REFERENCE.



## Before You Start

- ⚠ Please read all instructions carefully.
- ⚠ Retain instructions for future reference.
- ⚠ Separate and count all parts and hardware.
- ⚠ Read through each step carefully and follow the proper order.
- ⚠ We recommend that, where possible, all items are assembled near to the area in which they will be placed in use, to avoid moving the product unnecessarily once assembled.
- ⚠ Always place the product on a flat, steady and stable surface.
- ⚠ Keep all small parts and packaging materials for this product away from babies and children as they potentially pose a serious choking hazard.

**Warm tip: for your safety, please read the manual carefully before use!**

Please note: This product can not replace wheelchair, only for indoor use. It is strictly prohibited to carry out forced driving on uneven and dangerous roads, otherwise it may cause irreparable injury to the product and damage to the user.



## FOREWORD

### Dear users:

Before using the toilet chair, please read the operation manual carefully and use it correctly.

Please take good care of this manual after reading , so that you can consult it at any time.

## PRODUCT DESCRIPTION

This multifunctional chair is mainly made of aluminum alloy, with adjustable height caster, anti toppling foot pedal and anti-skid push hand, so that the elderly can use it safely. More humanized, it can solve the problems often encountered in the life of the elderly and pregnant women.

The weight-bearing is up to 100kg (220lbs), which greatly improves the dependence and demand of the elderly and pregnant women on the toilet chair, and makes your elderly and wife have a more comfortable home.

## PRODUCT SIZE


L \* W \* H: 82cm x 53cm x (87cm-92cm) / 32.3" x 21" x (34.3"-36.2")

Sitting height (including cushion): 48cm-53cm / 19"-21"


Net weight: 7.3kg / 16.2lbs




## INSTALLATION METHOD

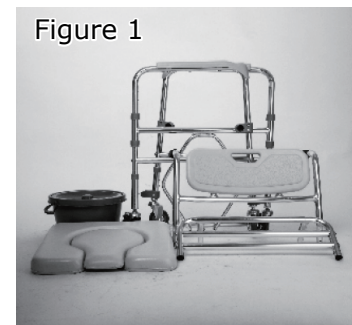


1. Open the box, take out the main frame and accessories, as shown in Figure 1.

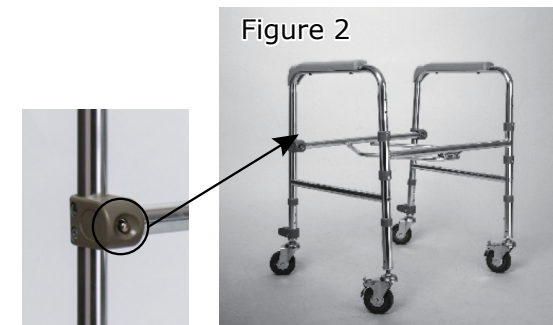

2. Install the rack, and the marbles shall be ejected. As shown in Figure 2.

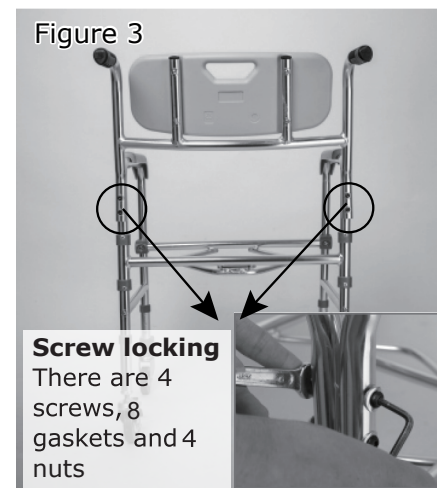


3. Install the backrest on the chair, align the hole position, and tighten the screw, as shown in Figure 3.

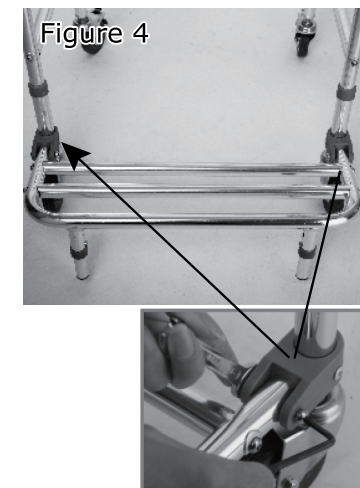

4. Install the foot pedal, align the hole position, and tighten the screw, as shown in Figure 4.

Figure 5

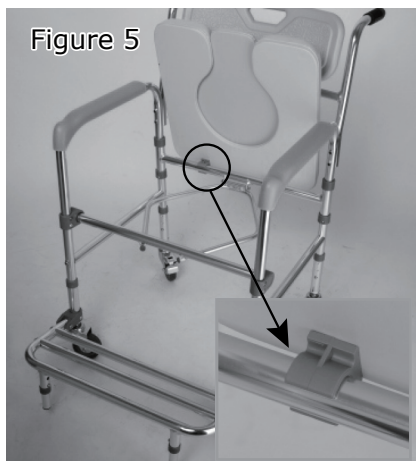

5. Put the cushion and clamp it. As shown in Figure 5.

Figure 6




6. Put on the toilet as shown in Figure 6.


Figure 7



7. Put down the cushion and you can use it. As shown in Figure 7.

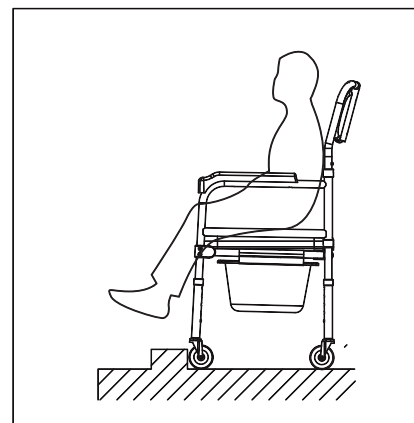

PROHIBIT

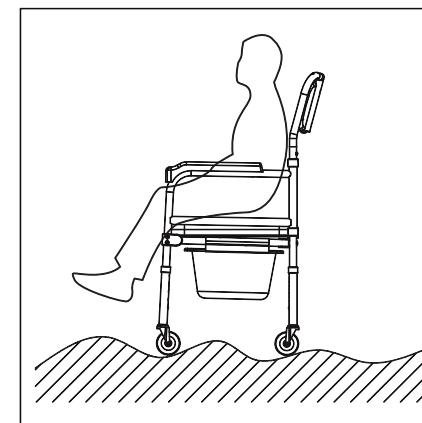

PROHIBIT

### Please note:

It is strictly prohibited to force the use of wheeled toilet chairs on uneven and bumpy roads or outdoors, otherwise it may cause irreparable injury to the product and damage to the user.