



COSTWAY®



COSTWAY
www.costway.com

USER'S MANUAL

Wheel Standard Walker with Arm Rest SP0602

THIS INSTRUCTION BOOKLET CONTAINS **IMPORTANT** SAFETY INFORMATION. PLEASE READ AND KEEP FOR FUTURE REFERENCE.

Please give us a chance to make it right and do better!

Contact our friendly customer service department for help first.

Replacements for missing or damaged parts will be shipped ASAP!

Follow Costway

Visit us: www.costway.com



Contact Us!

Do NOT return this item.

Contact our friendly customer service department for help first.



Before You Start

- ⚠ Please read all instructions carefully.
- ⚠ Retain instructions for future reference.
- ⚠ Separate and count all parts and hardware.
- ⚠ Read through each step carefully and follow the proper order.
- ⚠ We recommend that, where possible, all items are assembled near to the area in which they will be placed in use, to avoid moving the product unnecessarily once assembled.
- ⚠ Always place the product on a flat, steady and stable surface.
- ⚠ Keep all small parts and packaging materials for this product away from babies and children as they potentially pose a serious choking hazard.



Dear Customers:

- Thank you for using our standard walker.
- Before using it, please read the instruction manual carefully and store it in a safe place for future reference.
- If you need any further assistance, please contact us.


1 Specification

Height (adjustable): 30-36.6"/76-93cm
 Suitable for Human Height : 4'7"-6'1" / 140 - 185 cm
 Bearing Capacity: 440lb/200kg
 Folding Size: 27.16"x30"x4.7"/69cmx76cmx12cm
 Unfolding Size: 25.6"x20"x(30-36.6")/65cmx20cmx(76-93cm)
 Fixed Casters size: 5 inch
 Material: aluminum alloy
 Type: walking frame
 Applicable People: people with limited mobility, disabled people, rehabilitation workers, etc
 Advantages: easy to carry, foldable, 2 wheels are light and easy to move
 Applicable Occasions: home, senior apartment, nursing home, hospital
 Additional Features: height adjustable

2 Safety Precautions

Before use follow the precautions below:

1. Ensure all the legs on the walker are set to equal heights and that this is the correct height for the user.
2. Ensure the snap buttons have been securely snapped into the slot individually (The snap button is completely exposed outside the slot, and an audible "click" will be heard. Or the walker cannot remain steady.

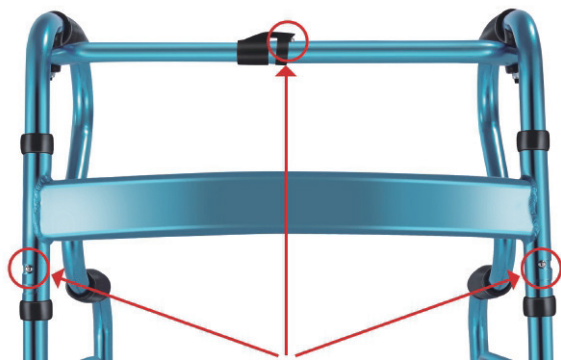

3. Do not exceed the weight limit of 440 lbs / 200 kg
4. Inspect the walker regularly for signs of damage or wear (see care and maintenance section)

3 Usage

3.1 Setting the walker



①



②



③




④



⑤

1. Take out the walker and all of the accessories from the carton.
2. Swing both side frames out away from the front cross brace until both snap buttons on them are fully exposed and stuck into the slot (If it cannot be stuck, turn the side frame back and forth until the hidden button on the side frame is completely exposed and stuck in the slot).
3. Pull the bar outward with both hands to fully expose the snap button on the beam and stuck it into the slot. Or the walker cannot remain steady
4. Connect the left side front wheeled tube with the walker, press down snap button and slide on, line up pin and hole so pin pops through walker hole completely. Repeat for the right side front wheeled tube.
5. Position the left side forearm platform up, connect it with the walker, press down snap button and slide on, line up pin and hole so pin pops through walker hole completely. In the meantime, fit the clip on the frame tube. Repeat for the right side forearm platform.
6. The installation is complete.

3.2 Folding the walker

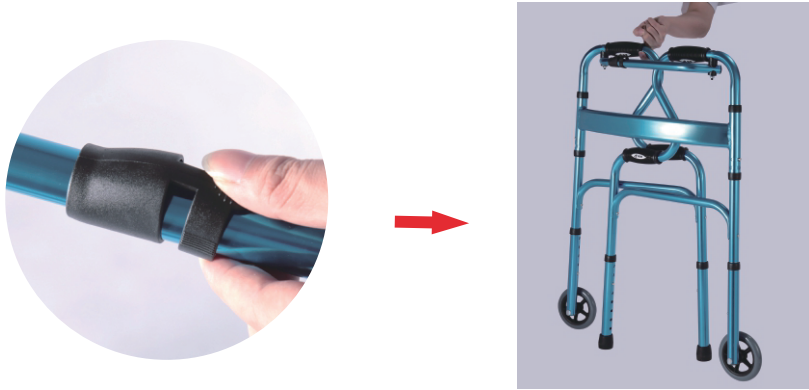

Find the plastic piece on the top of the front beam, there is a snap button under it. Press it down and fold the side frames in toward the front cross brace.

3.3 To adjust the height




1. Press the snap button, extend or retract the leg to raise or lower the height.
2. Adjust the leg extension to the desired height and allow the snap button to snap into the height adjustment hole. (An audible "click" will be heard)
3. Make sure that the walker is level.
Repeat for the other legs, ensuring all legs are securely in position and all adjusted to the same height.

3.4 Fixing / Interactive adjustment




Fixed Movement

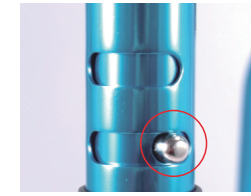

Interactive Movement

1. Fixed movement mode: press both the snap buttons into upper hole. Will be moved forward with lifting walker. **Fixed mode is recommended when using 2 wheels.**
 2. Interactive movement mode: press both the snap buttons into lower hole. Will be moved forward one side at a time, in step with the user. It allows the side frames of the walker to have a certain left and right movement space. **ONLY Suitable for using 4 walker sticks /legs, not for using wheels.**
- Direction for use: Move one side forward in step at a time. Lift the left side of the walker, move the left foot forward, and then lift the right side, and move the right foot forward.
3. When using the walker, please ensure that both left and right snap buttons are set to the same mode (on the same horizontal line), or the walker will not be balanced.

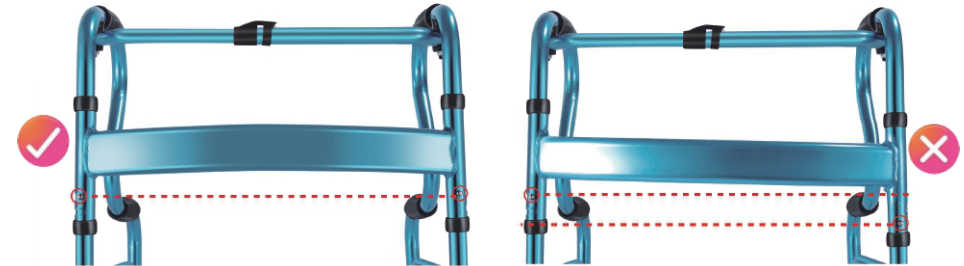


4 Care And Maintenance

Maintenance should be done at least every six months.
DO NOT use an abrasive detergent, cleaner or cloth to clean or dry the walker.

1. Regular cleaning of the walker is recommended. Wipe all surfaces with a damp cloth and allow all surfaces to dry before use.
2. Regular inspection of parts including hardware, brackets, and plastics, for deformation, corrosion, breakage, wear or compression is recommended. Replace the walker if any of these conditions exist.
3. If rubber tips are worn or cracked, replace the rubber tips IMMEDIATELY

5 How to Deal With the Unstable Problem With The Walker

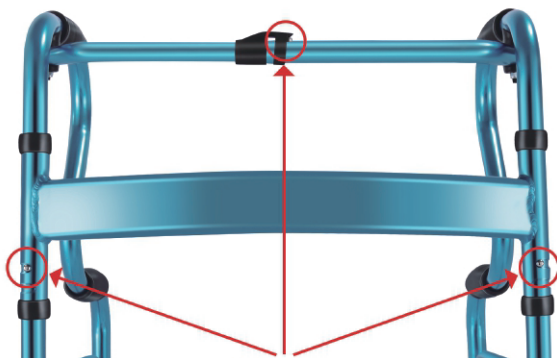

There are three snap buttons related to the stability of the walker. Please make sure that these three buttons are fully exposed and stuck into the slot.


No. 1

Pull the bar outward with both hands to fully expose the snap button on the beam and stuck it into the slot. Or the walker cannot remain steady.




No. 2

When using the walker, please ensure that both left and right snap buttons are set to the same mode (on the same horizontal line), or the walker will not be balanced.

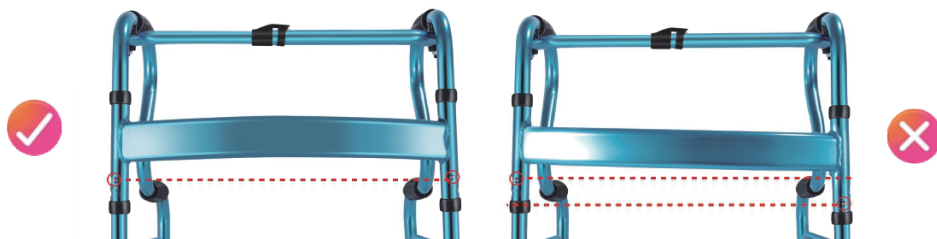

No. 3

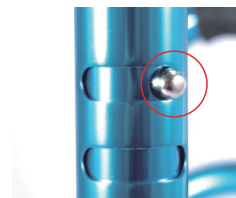

Fixed Movement

Fixed Movement Mode:
Press both the snap buttons into upper hole. Will be moved forward with lifting walker. *When using walker with 2 wheels, please use the fixed mode, otherwise the walker will be difficult to control.*



Interactive Movement

Interactive Movement Mode:
Press both the snap buttons into lower hole. Will be moved forward one side at a time, in step with the user. It allows the side frames of the walker to have a certain left and right movement space. *ONLY Suitable for using 4 walker sticks /legs, not for using wheels.*
Direction for use: Move one side forward in step at a time. Lift the left side of the walker, move the left foot forward, and then lift the right side, and move the right foot forward.

Welcome to visit our website and purchase our quality products!

With your inspiring rating, COSTWAY will be more consistent to offer you
EASY SHOPPING EXPERIENCE, GOOD PRODUCTS and EFFICIENT SERVICE!

US office: Fontana, California

UK office: Ipswich



- * Reward Points
- * Exclusive Customer Service
- * Personalized Recommendations
- * Permanent Shopping Cart
- * Order History

5 Star Rating



We would like to extend our heartfelt thanks to
all of our customers for taking time to assemble
this product and giving us valuable feedbacks.

